

TEAC

Are You Ready for DICOM Integration?



Medical Grade
Digital Video Recorder

UR-4MD

Designed for the mission critical OR environment

Full HD
1920 x 1080
Video and Photo recording

**Advanced
Reliability**
w/ Fail-safe recording


**DICOM
Gateway***
Modality Worklist support

3.5 inch
Built-in LCD Display

**Simultaneous recording
USB & iHDD**

Plug & Play
Automatic configuration

* DICOM Gateway is an optional upgrade.

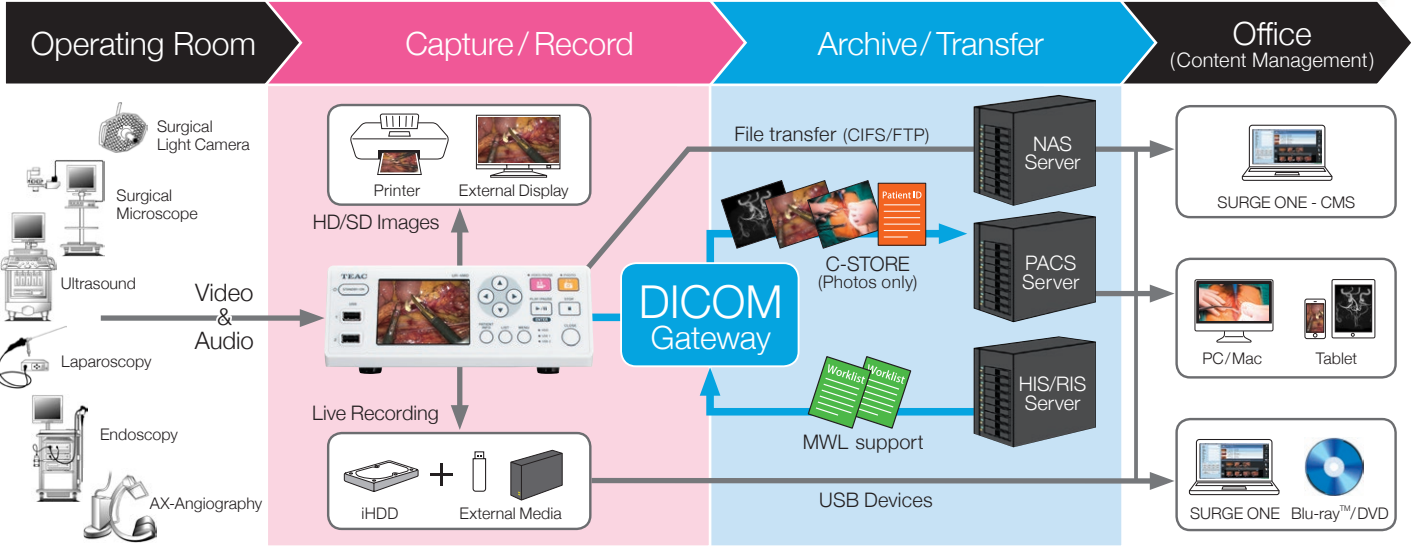

Key Features

- Crystal Clear, Full HD 1920 x 1080, Video & Photo Recording
- 3.5" built-in LCD Display for Live-view/Playback/Navigation
- Auto-Signal Detection & Auto-Resolution setting
- Safe Recording to iHDD, USB devices and Network drive simultaneously.
- Networking (FTP & CIFS) RJ45 10/100/1000 Ethernet
- Remote control via Scope, Foot-pedal, Barcode-reader or Keyboard (Optional)
- PACS connectivity: DICOM Storage & MWL support (Optional)
- Medical Standards

Specifications

Storage & Display			
Recordable Devices	USB Flash Drive, USB Hard Drive Internal Hard Drive (500GB) *1TB Optional Upgrade available		
Built-in LCD Display	3.5 inch Full Color		
Recording Format (FAT32)			
Photo	Compression: JPEG, DICOM; File: JPG, DCM		
Video	Compression: MPEG-4 AVC/H.264; File: MOV		
Input/Output			
Video	Input/Output	SDI (BNC) x1, DVI (Single Link) x1, S-Video x1, Composite (BNC) x1	3G/HD-SDI: 1080p/1080i/720p, SD-SDI: 480i/576i
	Output	SDI Thru x1	High-Definition, NTSC, PAL
	Input/Output	SDI (BNC) x1, Stereo mini jack x1	* SDI video output terminal used in common, LINE
	Input	Stereo mini jack x1	MIC
Other Interfaces	Input	RJ-45 (10BASE-T/100BASE-TX/1000BASE-T) x1	
	Output	USB 2.0 Type A x3	2 front connectors for Recordable Devices, 1 rear connector for Optional Accessories
	Output	USB 2.0 Type B x1	1 rear connector for accessing the internal Hard Drive (on UR-4MD)
General	Power supply	AC 100V - 240V, 50/60Hz	
Dimensions (WxDxH)	Approx. 8 1/4 x 9 1/4 x 3 3/8 in (210.5 x 235 x 88.5 mm)		
Weight	Approx. 5.7 lbs. (2.6 kg)		
Certifications	Safety	UL60601-2003, CE, IEC60601-1, CSA60601-1, EN60601-1	
Category	EMC	IEC60601-1-2, EN60601-1-2, ICES-003, FCC Class B	
What's in the box	FDA Class 1, ANVISA, RoHS, CCC		
Accessories (Optional)	TEAC UR-4MD, Power cords, Cable clamp, Owner's manuals		
	Remote control unit, Foot-switch, Keyboard, Barcode / Card Reader		

*Blu-ray™ Disc, "Blu-ray™" and "Blu-ray™" logo are trademarks of Blu-ray™ Disc Association. Other company names, product names and logos in this document are the trademarks or registered trademarks of their respective owners.
 * Features and specifications are subject to change without notice.

Product Overview

UR-4MD is TEAC's latest high-definition medical video recorder designed for quick set up and ease of use, yet includes a robust list of features and capabilities. With auto sensing video connections and automatic video resolution detection, set up is quick and easy. The front control panel allows users to intuitively record videos and photos in crystal clear HD. A built-in hard drive is included for reliable recording, data back-up and users can also simultaneously record to USB memory and USB hard drives to archive, share and play back on other devices and locations. UR-4MD is also Network capable. With an optional upgrade, DICOM integration is made easy.

Recording Times (for a 500GB Hard Drive)

High-Definition 1920 x 1080		
12 Mbps	18 Mbps	24 Mbps
90 hours	60 hours	45 hours
Standard-Definition 640 x 480		
5 Mbps	8 Mbps	10 Mbps
180 hours	135 hours	100 hours

Times are approximate.

UR-4MD Digital Video Recorder



TEAC America, Inc.,
 1834 Gage Road, Montebello, California, 90640, U.S.A.
 Tel.: +1-323-726-0303 Fax: +1-323-727-7688
 Email: mrsales@teac.com http://teac-ipd.com/